



## Forest Dental Component Catalogue

LIGHT SERVICE  
ITEMS **1**

HANDPIECE  
CONTROLS **2**

ARM  
SYSTEMS **3**

SYRINGES **4**

FOOT  
CONTROLS **5**

VACUUM **6**

HOLDERS **7**

VALVES **8**

ACCESSORIES **9**

FITTINGS **10**

TUBING **11**

SERVICE  
ITEMS **12**

PARTS  
BREAKDOWN **13**

### Phantom Head Dental

950 Manchester Road  
Linthwaite, Huddersfield, HD7 5QS

Tel/Ans: 01484 843 100  
Fax: 01484 843 888

Email: [sales@phantomhead.com](mailto:sales@phantomhead.com)  
Website: [www.phantomhead.com](http://www.phantomhead.com)



# Vacuum

6

## UNIVERSAL ASSISTANT KITS (325°F)



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- Autoclavable universal central vacuum lever valve with tubing, 5'
- Autoclavable saliva ejector tip valve with tubing, 5'
- Forestar-3 standard syringe with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-600	1400-500	1400-800	1400-900
HORIZONTAL SWING	35	1400-601	1400-501	1400-801	1400-901
VERTICAL SWING	35	1400-602	1400-502	1400-802	1400-902
w/o HOLDER BAR		1400-603	1400-503	1400-803	1400-903

## STANDARD ASSISTANT KIT (325°F)



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- Autoclavable standard central vacuum lever valve with tubing, 5'
- Autoclavable saliva ejector tip valve with tubing, 5'
- Forestar-3 standard syringe with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-604	1400-504	1400-804	1400-904
HORIZONTAL SWING	35	1400-605	1400-505	1400-805	1400-905
VERTICAL SWING	35	1400-606	1400-506	1400-806	1400-906
w/o HOLDER BAR		1400-607	1400-507	1400-807	1400-907



# Vacuum

6

## UNIVERSAL ASSISTANT KITS



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- White universal central vacuum lever valve with tubing, 5'
- Saliva ejector tip valve with tubing, 5'
- Forestar-3 standard syringe with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-608	1400-508	1400-808	1400-908
HORIZONTAL SWING	35	1400-609	1400-509	1400-809	1400-909
VERTICAL SWING	35	1400-610	1400-510	1400-810	1400-910
w/o HOLDER BAR		1400-611	1400-511	1400-811	1400-911

## STANDARD ASSISTANT KIT



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- White standard central vacuum lever valve with tubing, 5'
- Saliva ejector tip valve with tubing, 5'
- Forestar-3 standard syringe with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-612	1400-512	1400-812	1400-912
HORIZONTAL SWING	35	1400-613	1400-513	1400-813	1400-913
VERTICAL SWING	35	1400-614	1400-514	1400-814	1400-914
w/o HOLDER BAR		1400-615	1400-515	1400-815	1400-915



# Vacuum

6

## AUTOCLAVABLE UNIVERSAL CENTRAL VACUUM KITS (325°F)



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- Autoclavable central vacuum lever valves with tubing, 5'
- Autoclavable saliva ejector tip valve with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-278	1400-078	1400-878	1400-978
HORIZONTAL SWING	35	1400-276	1400-177	1400-877	1400-977
VERTICAL SWING	35	1400-267	1400-167	1400-876	1400-976
w/o HOLDER BAR		1400-274	1400-074	1400-875	1400-975

## AUTOCLAVABLE STANDARD CENTRAL VACUUM KIT (325°F)



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- Autoclavable central vacuum lever valves with tubing, 5'
- Autoclavable saliva ejector tip valve with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-277	1400-077	1400-874	1400-974
HORIZONTAL SWING	36	1400-263	1400-163	1400-873	1400-973
VERTICAL SWING	36	1400-264	1400-164	1400-872	1400-972
w/o HOLDER BAR		1400-282	1400-082	1400-871	1400-971



# Vacuum

6

## UNIVERSAL CENTRAL VACUUM KIT

### UNIVERSAL LEVER VALVE KIT



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- White universal central vacuum lever valves with tubing, 5'
- Saliva ejector tip valve with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-258	1400-158	1400-869	1400-969
HORIZONTAL SWING	35	1400-275	1400-175	1400-868	1400-968
VERTICAL SWING	35	1400-265	1400-165	1400-867	1400-967
w/o HOLDER BAR		1400-261	1400-161	1400-866	1400-966

## STANDARD CENTRAL VACUUM KIT



Used with central vacuum.

**INCLUDES**

- Central vacuum canister with disposable solids collector
- White standard central vacuum lever valves with tubing, 5'
- Saliva ejector tip valve with tubing, 5'
- Three position holder bar
- Mounting hardware



DESCRIPTION	MOUNT PAGE #	PART NO. STERLING	PART NO. GRAY	PART NO. BLACK	PART NO. SHADOW
FIXED POSITION	37	1400-257	1400-157	1400-865	1400-965
HORIZONTAL SWING	36	1400-253	1400-153	1400-864	1400-964
VERTICAL SWING	36	1400-254	1400-154	1400-863	1400-963
w/o HOLDER BAR		1400-260	1400-160	1400-862	1400-962



# Vacuum

6

## AUTOCLAVABLE VACUUM LEVER VALVES (325°F)



STANDARD



Autoclavable valve features finger tip on/off and flow control with unique lever operated design. Swivel quick disconnect allows 360 degrees rotation of entire valve and adapts to any 1/2" I.D. central vacuum tubing. After valve is installed in tubing, it can be snapped apart at the swivel to allow autoclaving of valve assembly (lever, body, and o-rings). Universal valve has one piece adapter added for use with universal type handpiece holders. Both valves accept 7/16" O.D. tips.




UNIVERSAL

DESCRIPTION	PART NO.
STANDARD	1100-250
UNIVERSAL	1100-251

**Note:** Replacement o-rings 0015-313 (#13), see pg. 78

## STANDARD VACUUM LEVER VALVES

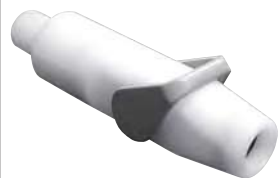

STANDARD



Features finger tip on/off and flow control with unique lever operated design. Swivel quick disconnect allows 360 degrees rotation of entire valve and adapts to any 1/2" I.D. central vacuum tubing. Universal valve has one piece, clear anodized satin finish adapter added for use with universal type handpiece holders. Both valves accept 7/16" O.D. tips.



UNIVERSAL

DESCRIPTION	PART NO.
STANDARD (WHITE)	1100-230-WHT
STANDARD (BLACK)	1100-230-BLK
UNIVERSAL (WHITE)	1100-231-WHT
UNIVERSAL (BLACK)	1100-231-BLK

**Note:** Replacement o-ring 0015-013 (#13), see pg. 79

## VACUUM ROTARY VALVES



STANDARD



Features on/off and flow control with rotary type valve. Swivel quick disconnect allows 360 degrees rotation of entire valve and adapts to any 1/2" I.D. central vacuum tubing. Universal valve has one piece, clear anodized satin finish adapter added for use with universal type handpiece holders. Both valves accept 7/16" O.D. tips.



UNIVERSAL

DESCRIPTION	PART NO.
STANDARD (BLACK)	1100-106
UNIVERSAL (BLACK)	1100-111

**Note:** Replacement o-ring 0015-013 (#13), see pg. 79



# Vacuum

6

## AUTOCLAVABLE SALIVA EJECTORS (325°F)



SALIVA EJECTOR VALVE



Features finger tip on/off and flow control with unique lever operated design. Swivel quick disconnect works with universal holders allowing 360 degrees rotation of entire valve and adapts to any 3/8" O.D. saliva ejector tubing. After valve is installed in tubing, it can be snapped apart at the base to allow autoclaving of valve assembly (lever, body and o-rings).

**DESCRIPTION**  
SALIVA EJECTOR VALVE  
SALIVA EJECTOR BODY

**PART NO.**  
1100-307  
1100-304



SALIVA EJECTOR BODY

**Note:** Replacement o-rings 0015-010 & 0015-308 (#8), see pg. 78 & 79

## SALIVA EJECTOR TIP VALVE ASSEMBLY



PUSH-ON TIP



Features rotary type on/off valve, with either push-on or screw-on style tip.



SCREW-ON TIP

DESCRIPTION	PART NO.
BLACK BODY WITH PUSH-ON TIP (BLACK)	1100-107
BLACK BODY WITH SCREW-ON TIP (BLACK)	1100-262
BLACK BODY WITH PUSH-ON TIP (GRAY)	1100-267
WHITE BODY WITH PUSH-ON TIP (GRAY)	1100-267-WHT

## SALIVA EJECTOR TIP ASSEMBLY



PUSH-ON TIP



Features either push-on or screw-on type tips.

**Note:** Available only in black.



SCREW-ON TIP

DESCRIPTION	PART NO.
WITH PUSH-ON TIP	1100-115
WITH SCREW-ON TIP	1100-261

## SALIVA EJECTOR ADAPTER



Adapts saliva ejector tip for use in central vacuum valve.

**DESCRIPTION**  
SALIVA EJECTOR ADAPTER

**PART NO.**  
1100-264



# Vacuum

6

## REPLACEMENT PARTS



**AUTOCLAVABLE SALIVA EJECTOR TIP (325°F)**  
For use with saliva ejector assemblies.

**DESCRIPTION**  
AUTOCLAVABLE (GRAY)  
AUTOCLAVABLE (BLACK)

**PART NO.**  
1100-308  
1100-108

## SALIVA EJECTOR TIP



For use with saliva ejector assemblies.

**DESCRIPTION**  
SCREW-ON (BLACK)

**PART NO.**  
1100-258

## SALIVA EJECTOR SCREEN



Fits push-on or screw-on style saliva ejector assemblies.

**DESCRIPTION**  
PLASTIC SCREEN

**PART NO.**  
0010-044

## BASE



Saliva ejector base adapts to any 3/8" O.D. saliva ejector tubing. The central vacuum base adapts to any 1/2" I.D. central vacuum tubing. Both bases allow for 360 degrees rotation of valve.

SALIVA EJECTOR  
CENTRAL VACUUM

**DESCRIPTION**  
SALIVA EJECTOR  
CENTRAL VACUUM

**PART NO.**  
1100-314  
1100-320

## LEVER



Allows for on/off and flow control with use of finger tip. Both levers are autoclavable (325°F).

**DESCRIPTION**  
SALIVA EJECTOR (GRAY)  
SALIVA EJECTOR (BLACK)  
CENTRAL VACUUM (GRAY)  
CENTRAL VACUUM (BLACK)

**PART NO.**  
1100-311  
1100-311-BLK  
1100-248  
1100-247

## PLUG



Stops suction and allows saliva ejector and central vacuum tubings to be placed in universal holder after valve has been removed from quick disconnect base.

**DESCRIPTION**  
PLUG

**PART NO.**  
1400-136



# Vacuum

6

## CENTRAL VACUUM CANISTERS

STANDARD VACUUM CANISTER, 1/2" DRAIN



Features removable solids collector, two high velocity ports for 1/2" I.D. tubing, one saliva ejector port for 3/8" I.D. tubing, optional water drain port for 1/8" I.D. tubing, extra port plugs, mounting brackets and hardware.

DESCRIPTION	PART NO.
CANISTER ASSEMBLY (GRAY)	1400-111
CANISTER ASSEMBLY (BLACK)	1400-011
CAP & O-RING ONLY (GRAY)	1400-155
CAP & O-RING ONLY (BLACK)	1400-055
O-RING ONLY	0015-033
SOLIDS COLLECTOR (GREEN) (Pkg. of 10)	0020-100
CV HORIZONTAL BRACKET (GRAY)	1400-128
CV HORIZONTAL BRACKET (BLACK)	1400-128-BLK
CV VERTICAL BRACKET (GRAY)	1400-129
CV VERTICAL BRACKET (BLACK)	1400-129-BLK

## VACUUM CANISTER, 1/2" DRAIN WITH LOCKING LID



Features same as above except with bayonet type locking lid.

DESCRIPTION	PART NO.
CANISTER ASSEMBLY (GRAY)	1400-190
CANISTER ASSEMBLY (BLACK)	1400-090
CAP & O-RING ONLY (GRAY)	1400-193
CAP & O-RING ONLY (BLACK)	1400-093
O-RING ONLY	0015-033
SOLIDS COLLECTOR (GREEN) (Pkg. of 10)	0020-100
CV HORIZONTAL BRACKET (GRAY)	1400-128
CV HORIZONTAL BRACKET (BLACK)	1400-128-BLK
CV VERTICAL BRACKET (GRAY)	1400-129
CV VERTICAL BRACKET (BLACK)	1400-129-BLK

## VACUUM CANISTER, 1/2" REAR DRAIN WITH LOCKING LID



Features removable solids collector, two high velocity ports for 1/2" I.D. tubing, one saliva ejector port for 3/8" I.D. tubing, extra port plugs, bayonet type locking lid and rear drain. For use with Vacuum Telescoping Arm System.

DESCRIPTION	PART NO.
CANISTER ASSEMBLY (WHITE)	1400-133
CANISTER ASSEMBLY (BLACK)	1400-133-BLK
SOLIDS COLLECTOR (BLUE) (Pkg. of 10)	0020-101

## CENTRAL VACUUM CANISTER PLUGS



DESCRIPTION	PART NO.
SALIVA EJECTOR VALVE PORT (WHITE)	1400-008
CENTRAL VACUUM VALVE PORT (WHITE)	1400-017
CENTRAL VACUUM VALVE PORT (BLACK)	1400-017-BLK



## Vacuum

6

### SOLIDS COLLECTOR



For use with central vacuum canisters. Green fits 1400-111 standard and 1400-190 locking cannister assemblies. Blue fits 1400-133 5/8" rear drain cannister assembly.

**DESCRIPTION**  
 PACKAGE OF 10 (GREEN)  
 PACKAGE OF 10 (BLUE)

**PART NO.**  
 0020-100  
 0020-101

### SOLIDS COLLECTOR W/ LID



For use with central vacuum canisters.  
 Small - fits 1400-111 standard and 1400-190 locking cannister assembly.  
 Large - fits 1400-133, 5/8" rear drain cannister assembly.

**DESCRIPTION**  
 EACH (CLEAR: SMALL)  
 EACH (CLEAR: LARGE)

**PART NO.**  
 1400-191  
 1400-192

### POST MOUNT ADAPTER



Allows either central vacuum canister, 1400-111 or 1400-190, to be mounted to a 2" O.D. post.

**DESCRIPTION**  
 POST MOUNT ADAPTER

**PART NO.**  
 0023-436

### SALIVA EJECTOR ADAPTERS

ADAPTER, 1/4" I.D. X 7/16" I.D.



Adapts 3/8" O.D. saliva ejector tubing to vacuum valve using 7/16" O.D. tips.

**DESCRIPTION**  
 SALIVA EJECTOR ADAPTER

**PART NO.**  
 0016-114

## Phantom Head Dental

950 Manchester Road  
 Linthwaite, Huddersfield, HD7 5QS

Tel/Ans: 01484 843 100  
 Fax: 01484 843 888

Email: sales@phantomhead.com  
 Website: www.phantomhead.com



## Vacuum

6

### SALIVA EJECTOR ADAPTERS

ADAPTER, 1/4" I.D. X 1/4" I.D.

Adapts 3/8" O.D. saliva ejector tubing to saliva ejector tip.



**DESCRIPTION**  
SALIVA EJECTOR ADAPTER

**PART NO.**  
0016-214

### IN-LINE MANUAL ON/OFF VALVE



Manual operated vacuum on/off valve may be added in-line to 1/2" I.D. central vacuum tubing. Enables easy removal of vacuum canister lid without shutting vacuum down.

**DESCRIPTION**  
MANUAL ON/OFF VALVE

**PART NO.**  
1100-244

### ORAL DRY CUP



Autoclavable polysulfone (beige). Fits all vacuum valves with 7/16" dia. tips.

**DESCRIPTION**  
ORAL DRY CUP

**PART NO.**  
1400-010

### STANDARD HVE DISPOSABLE TIP



Made of lightweight polyethylene, 7/16" O.D. x 6" long. Packaged in quantities of 100.

**DESCRIPTION**  
PACKAGE OF 100

**PART NO.**  
1400-020

## Phantom Head Dental

950 Manchester Road  
Linthwaite, Huddersfield, HD7 5QS

Tel/Ans: 01484 843 100  
Fax: 01484 843 888

Email: sales@phantomhead.com  
Website: www.phantomhead.com



# Vacuum

6

## CUSPIDOR STRAINER



Molded white plastic. Fits Forest cuspidor bowl.

**DESCRIPTION**  
CUSPIDOR STRAINER

**PART NO.**  
0010-129

## VACUUM DRAIN CUSPIDOR KIT



May be used to adapt all models of gravity drain cuspidor to central vacuum system.

**INCLUDES**

- Vacuum drain valve
- Flush nipple
- Filter assembly
- Tubing
- Installation hardware

**DESCRIPTION**  
FOR 5/8" I.D. DRAIN  
FOR 7/8" I.D. DRAIN

**PART NO.**  
1400-116  
1400-117

## AIR/WATER SALIVA EJECTOR



Air/water operated saliva ejector.

**INCLUDES**

- Vacuum generator
- Saliva ejector tip with tubing, 6'
- Toggle valve
- 1/2" I.D. drain tubing, 9'
- Installation hardware

**DESCRIPTION**  
STERLING  
GRAY  
BLACK  
SHADOW

**PART NO.**  
1400-285  
1400-085  
1400-885  
1400-985

## SALIVA EJECTOR HOSE ASSEMBLY



Saliva ejector assembly adapts to any vacuum valve using 7/16" O.D. tips.

**INCLUDES**

- Saliva ejector tip
- 7/16" O.D. adapter
- Tubing, 3'

**DESCRIPTION**  
STERLING  
GRAY  
BLACK  
SHADOW

**PART NO.**  
1400-252  
1400-052  
1400-852  
1400-952



# Vacuum

6

## AUTOMATIC DUAL AIR VACUUM SYSTEM



Air-operated oral evacuation system activated by automatic holders.

**INCLUDES**

- White standard vacuum lever valve with tubing, 5'
- White saliva ejector tip valve with tubing, 5'
- 2 automatic holders
- Dual port venturi
- One quart (1 liters) waste container
- Air pilot valve



**DESCRIPTION**

STERLING  
GRAY  
BLACK  
SHADOW

**PART NO.**

1400-173  
1400-174  
1400-870  
1400-970

## MANUAL DUAL AIR VACUUM SYSTEM



Air-operated oral evacuation system activated by manual toggle valve.

**INCLUDES**

- White standard vacuum lever valve with tubing, 5'
- White saliva ejector tip valve with tubing, 5'
- 2 standard holders
- Dual port venturi
- One quart (1 liters) waste container
- 3 way toggle valve



**DESCRIPTION**

STERLING  
GRAY  
BLACK  
SHADOW

**PART NO.**

1400-183  
1400-184  
1400-861  
1400-961

